



GMM General Medical Merate S.p.A
www.gmmspa.com

Infotec



POLISTAT M

Universal radiographic stand

MEAN FEATURES

- Wall fixed stand with vertical slide device
- Vertically sliding arm stand, rotating on a central pivot supporting the radiation source assembly and the bucky
- Wide range of focal distance
- Manual movements with mechanical or electromagnetic braking

POLISTAT M is a wall fixed stand with a rotating arm supporting the radiation tube assembly and the bucky.

The unit allows easy exams performing with the patient standing, sitting or lying down.

Furthermore, it is especially suitable for basic radiological investigation and first aid, may be coupled with a radio transparent stretcher top (optional)

MAIN APPLICATIONS

- Emergency
- General diagnostic
- Skeletal and thoracic diagnostic








TECHNICAL DATA

STAND	
Column height	216 cm
ROTATING ARM	
Vertical travel	110 cm
Focus-floor minimum distance in horiz. position	50 cm
Movement type	manual
Balancing	counterweights
Rotation on the horizontal axis	+/- 120°
Braking device	electromechanical
ARM TUBE SUPPORT	
Manual rotation around the horizontal axis	-90/+180°
Angle position mechanically detected	0°
Incident angle indicator	goniometer
BUCKY SUPPORT	
Manual rotation on its axis	+/- 45°
Angular position mechanical detected	0°
VARIABLE FOCUS	
FFD range	98 ÷ 188 cm
Movement type	manual
FFD value visualization	graduate scale
CONTROLS	
Hand grips for movements	covered in synthetic rubber
Pushbuttons for brakes unlocking	in ergonomic position
Supplementary pushbutton for vertical movement brake	On bucky side
Power requirements	
Line voltage / frequency	230 (220/240) VAC 50/60Hz
Max absorbed power	300 VA
Weights	
Stand	380 kg

COMPLIMENTARY PARTS

	Manual collimator
	Potter-bucky
	Grid FFD=150 cm R=10:1 90 l/ inch
	Fixed height stretcher with radio transparent tabletop
	Variable height and tilting stretcher with radio transparent table

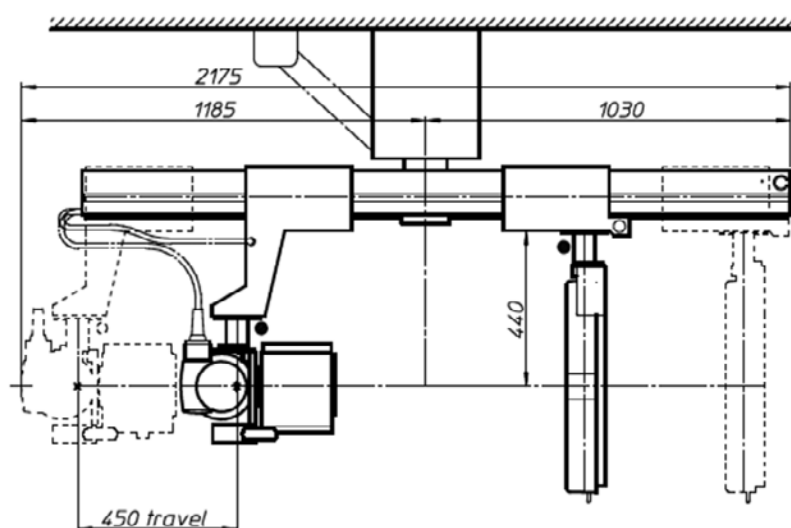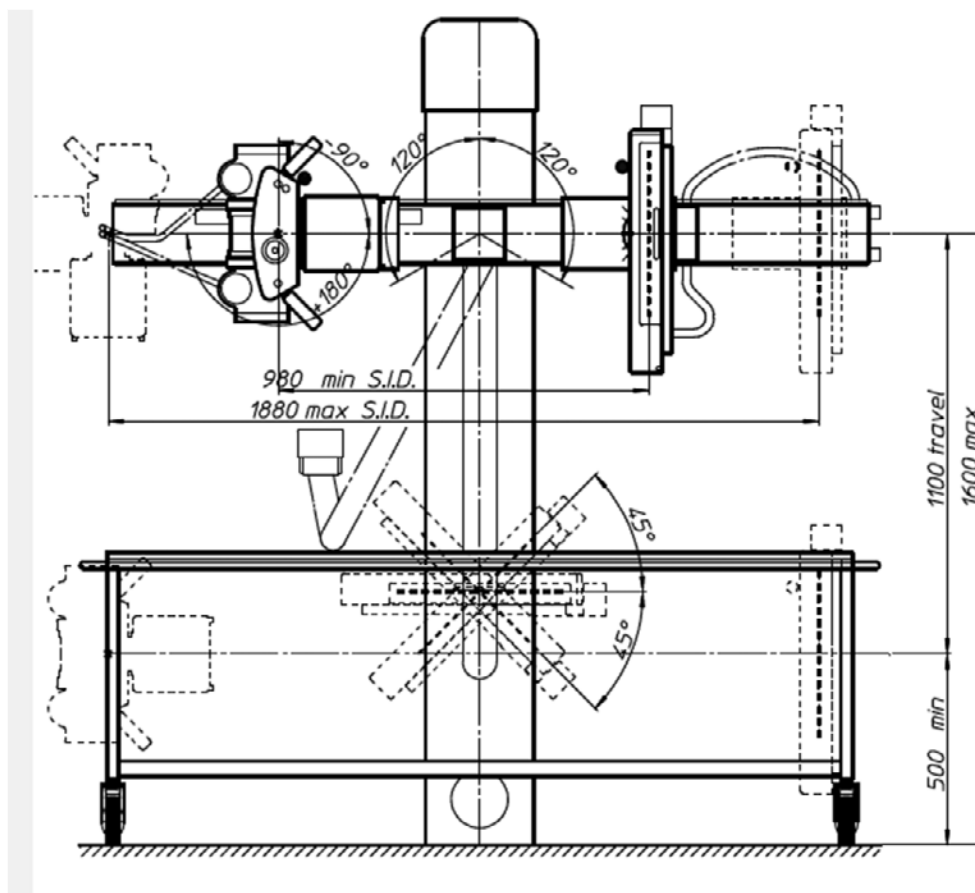

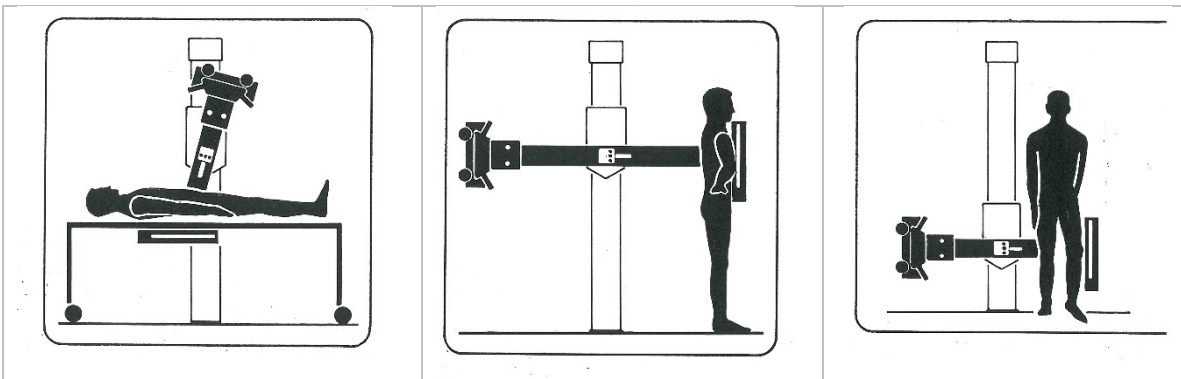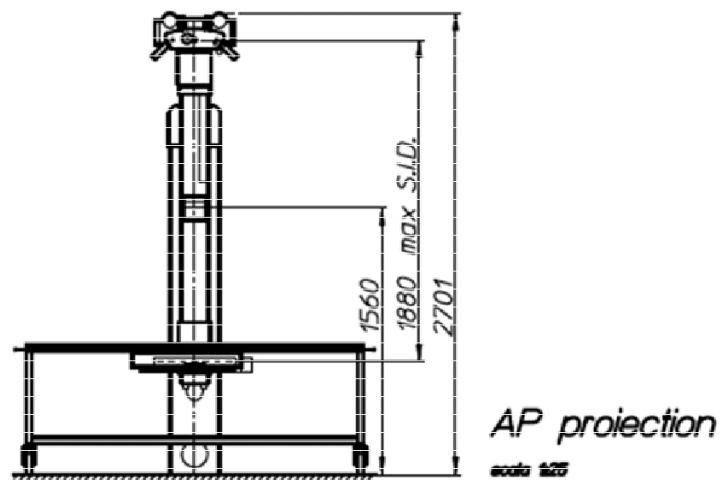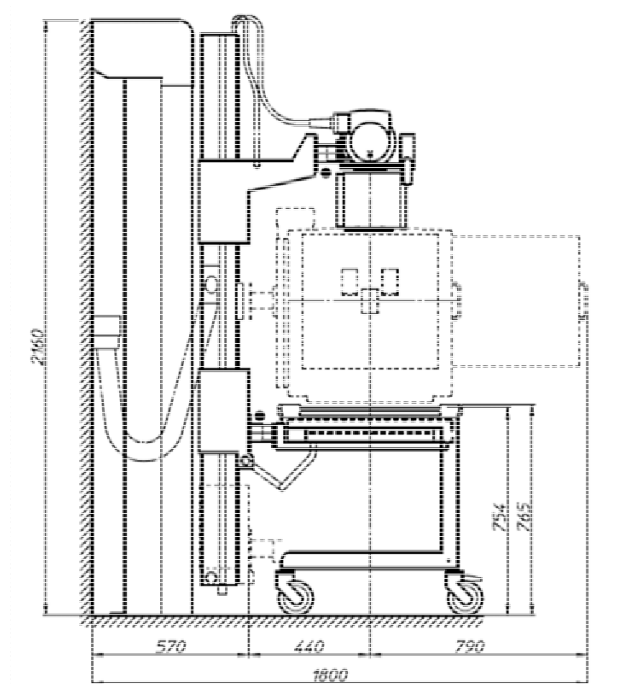



GENERAL MEDICAL MERATE S.p.A.

Via Partigiani, 25 - 24068 Seriate (BG) ITALIA

tel (0039) 035 4525311 - fax (0039) 035 297787

http: www.gmmspa.com e-mail: info@gmmspa.com

Conforme ai requisiti essenziali della DIRETTIVA DISPOSITIVI MEDICI 93/42CEE e sue integrazioni
In compliance with the essential requirements of MEDICAL DEVICES DIRECTIVE 93/42EEC and supplements